

M•Drain™ Single Step Drainage Catheter

The **M•Drain™ Single Step Drainage Catheter** is made of soft flexible polyurethane.



The catheter can be used for nephrostomy, abscess and other multipurpose drainage applications. Catheter placement is performed using Single Step or Seldinger technique. The easy-to-use, reliable safety string lock system provides strong fixation of the pigtail while maintaining patient comfort.

Hydrophilic coating, when activated, becomes very lubricous and reduces greatly the coefficient of friction. The decreased friction of the distal catheter tipping surface facilitates catheter introduction.

Large “skived” oval drainage holes that increase drainage capability.

Features:

- **Large “skived” oval side holes** promote increased flow and reduce clogging
- **Slim line reinforced tapered tip** with “Anti-Accordion” design for smooth entry
- **Unique leak proof pigtail locking system** to avoid mess and unnecessary catheter exchanges
- **Internal diameter of hub and catheter match** for maximum flow
- **Centimeter markings** to assist in accurate placement
- **Echogenic metal and plastic stiffeners** provide flexibility in catheter advancement
- **User friendly snap off locking system** provides improved patient comfort and catheter securement
- **Echogenic needle tip** with trocar stylet for easy, trackable insertion
- **Hydrophilic coating** reduces friction co-efficient facilitating catheter introduction
- **Easy-to-remove pigtail straightener** improves ease of use

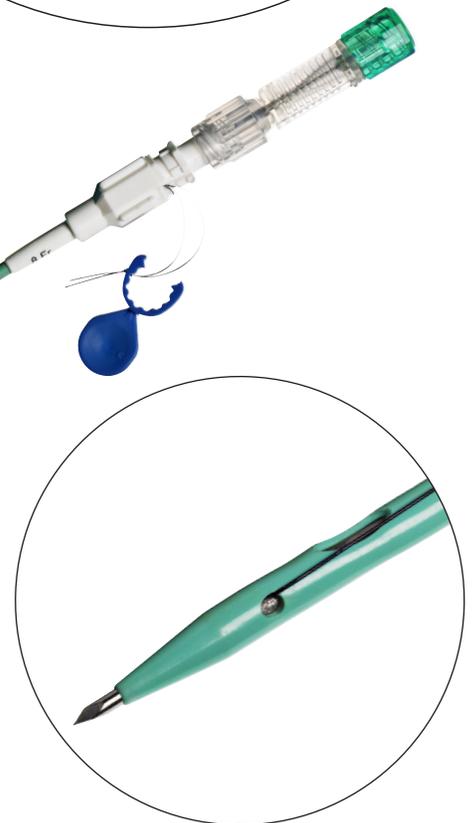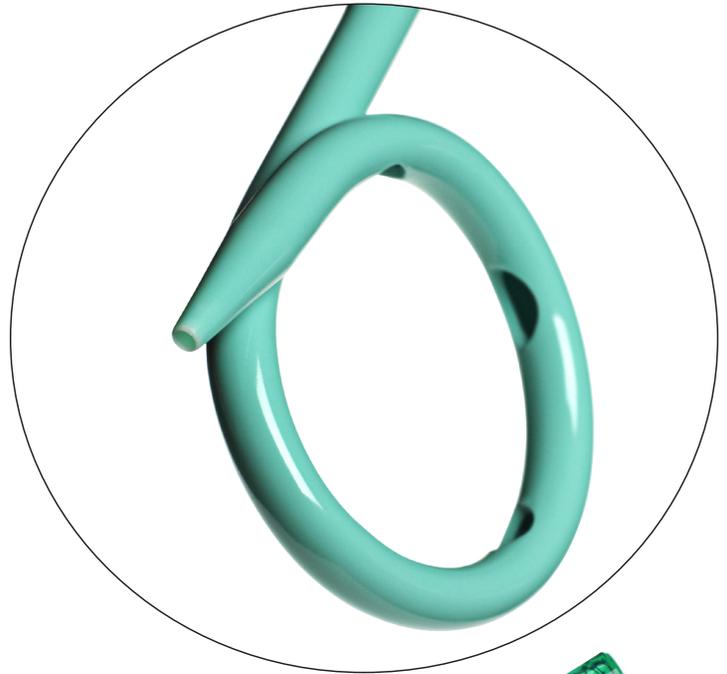

Full Matrix on page 2

M°Drain™ Single Step Drainage Catheter



All below numbers ending with "S" include flexible stiffener

M°Drain™ Catheter with locking pigtail

| Item no. | Size | Qty/box |
|-------------|-------------|---------|
| D11006020 | 6F x 20 CM | 5 |
| D11006025 | 6F x 25 CM | 5 |
| D11006035* | 6F x 35 CM | 5 |
| D11007015S* | 7F x 15 CM | 5 |
| D11007020S | 7F x 20 CM | 5 |
| D11007025S | 7F x 25 CM | 5 |
| D11007035S | 7F x 35 CM | 5 |
| D11008020S | 8F x 20 CM | 5 |
| D11008025S | 8F x 25 CM | 5 |
| D11008030S* | 8F x 30 CM | 5 |
| D11008035S | 8F x 35 CM | 5 |
| D11010015S* | 10F x 15 CM | 5 |
| D11010020S | 10F x 20 CM | 5 |
| D11010025S | 10F x 25 CM | 5 |
| D11010030S* | 10F x 30 CM | 5 |
| D11010035S | 10F x 35 CM | 5 |
| D11012025S | 12F x 25 CM | 5 |
| D11012035S | 12F x 35 CM | 5 |
| D11014025S | 14F x 25 CM | 5 |
| D11014035S | 14F x 35 CM | 5 |
| D11016025S | 16F x 25 CM | 5 |
| D11016035S | 16F x 35 CM | 5 |

M°Drain™ Catheter without locking pigtail

| Item no. | Size | Qty/box |
|-----------|-------------|---------|
| D10505015 | 5F x 15 CM | 5 |
| D10505020 | 5F x 20 CM | 5 |
| D10506015 | 6F x 15 CM | 5 |
| D10506020 | 6F x 20 CM | 5 |
| D10506025 | 6F x 25 CM | 5 |
| D10507020 | 7F x 20 CM | 5 |
| D10507025 | 7F x 25 CM | 5 |
| D10508020 | 8F x 20 CM | 5 |
| D10508025 | 8F x 25 CM | 5 |
| D10510020 | 10F x 20 CM | 5 |
| D10510025 | 10F x 25 CM | 5 |
| D10512025 | 12F x 25 CM | 5 |
| D10514025 | 14F x 25 CM | 5 |
| D10516025 | 16F x 25 CM | 5 |

* Special order product.

Please contact Customer Care for lead times.