

Sengstaken - Blakemore tubes

For emergency treatment
of bleeding oesophageal
varices

● A complete range of four lumen oesophageal tamponades for adults, children and paediatrics.


- One lumen for inflation of the gastric balloon.
- One lumen for inflation of the oesophageal balloon.
- One lumen for aspiration of secretions and blood, above the oesophageal balloon.
- One lumen for aspiration of gastric content.

● Made from rubber, with latex oesophageal and gastric sealing balloons.

● Radio - opaque marking for confirmation of tube placement.

● Sterile by Gamma Radiation

● 1 unit per box.

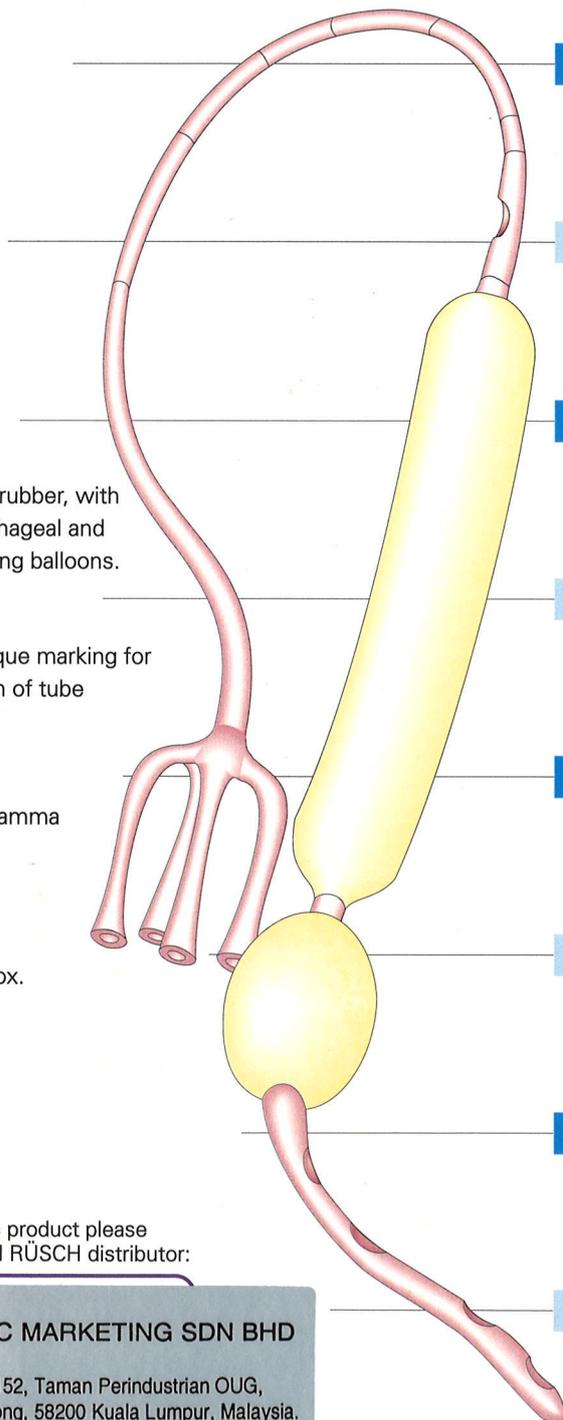

Product code	Product information
9160 14	Sengstaken four lumen tube 14Ch. Paediatric size Length: 960mm Oesophageal balloon length: 95mm Gastric balloon volume range: 50 - 500ml
9161 14	Sengstaken four lumen tube 14Ch. Adult size Length: 960mm Oesophageal balloon length: 185mm Gastric balloon volume range: 100 - 750ml
9162 20	Sengstaken four lumen tube 20Ch. Adult size Length: 960mm Oesophageal balloon length: 185mm Gastric balloon volume range: 100 - 750ml
9163 20	Sengstaken four lumen tube 20Ch. Adolescent size Length: 960mm Oesophageal balloon length: 115mm Gastric balloon volume range: 100 - 750ml
9164 16	Sengstaken four lumen tube 16Ch. Paediatric size Length: 960mm Oesophageal balloon length: 95mm Gastric balloon volume range: 50 - 100ml
9165 20	Sengstaken four lumen tube 20Ch. Adult size Gastric balloon only - 100 - 750ml Length: 960mm
9166 16	Sengstaken four lumen tube 16Ch. Adult size Length: 960mm Oesophageal balloon length: 185mm Gastric balloon volume range: 100 - 750ml
9167 18	Sengstaken four lumen tube 18Ch. Adolescent size Length: 960mm Oesophageal balloon length: 115mm Gastric balloon volume range: 100 - 750ml
9168 18	Sengstaken four lumen tube 18Ch. Adult size Length: 960mm Oesophageal balloon length: 185mm Gastric balloon volume range: 100 - 750ml

For supply of this product please contact your local RÜSCH distributor:

km KIMMAC MARKETING SDN BHD
(673056-V)

No. 5, Jalan 7/152, Taman Perindustrian OUG,
Batu 6, Jalan Puchong, 58200 Kuala Lumpur, Malaysia.
Tel: 03-7785 9293, 7785 9499, 7785 9255, 7785 9225
Fax: 03-7785 9252 E-mail: kimmac@po.jaring.my